


EP-TRACER EP RECORDING SYSTEM



Portable System

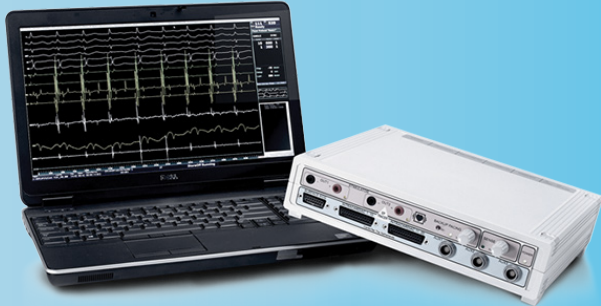

Fixed Installation



Mobile Cart

**Recorder
and Stimulator:**
2-in-1 system

**Excellent
Signal
Quality**

**Flexible
Configurations**

**Optimised
for Cryoablation**

REDEFINE YOUR EP LAB TREAT MORE PATIENTS

EP-TRACER is the electrophysiological measurement system with a fully-integrated stimulator.

The EP-TRACER Solution:


- 1 Recorder and Stimulator: 2-in-1 system**
 - Recorder with built-in stimulator controlled through the same software
 - All components fully-integrated
 - Critical information visualised on one screen
 - Less wiring/cable connection problems due to integrated design
- 2 Excellent Signal Quality**
 - Built-in stimulator helps reduce influence of external systems and provides optimal signal quality
- 3 Flexible Configurations**
 - Configurations include Mobile Cart, Portable System, Fixed Installation
 - Small footprint
 - Can be rapidly moved between laboratories and hospitals
 - Can turn any cath lab into fully functional EP lab
- 4 Optimised for Cryoablation**
 - Temperature curve display
 - Time to isolation
 - Visualisation of a CMAP trace to display phrenic nerve activity

EP-TRACER MODELS

- 3 amplifier models with different intracardiac channels (20, 52 or 84 ICs)
- All amplifier models allow connection of 12 surface ECGs and 6 auxiliary channels
- All models compatible with all system configurations (Portable System, Mobile Cart, Fixed Installation)
- Compatible with Cryoablation system, 3D mapping solutions

CORE PRODUCT

Amplifier with built-in stimulator

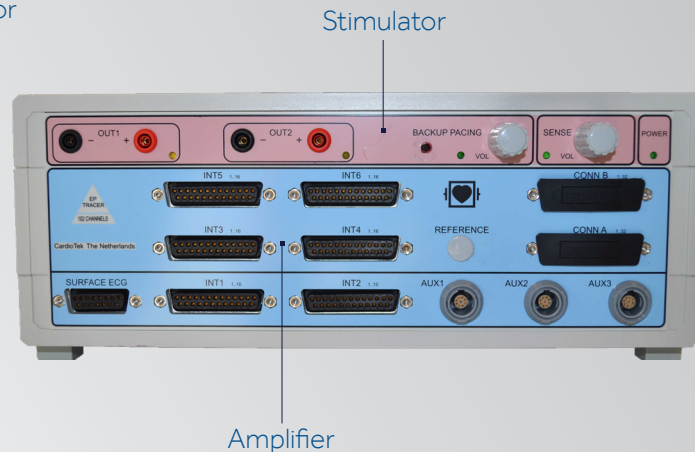

Brief Statement

See the device manual for detailed information regarding the instructions for use, indications, contraindications, warnings, precautions, and potential adverse events. For further information, contact your local Medtronic representative.

EP-TRACER is an electrophysiology measurement system used to acquire, filter, digitize, amplify, display, and record signals obtained during electrophysiological studies and related procedures.

The EP-TRACER received CE-Mark and FDA 510(k) clearance. The system is commercially available in Europe, Central and South America, Asia Pacific, Africa and the Middle East.

EP-TRACER EP recording system manufactured by CARDIOTEK BV, distributed by MEDTRONIC.

Medtronic

Medtronic International Trading Sarl
Route du Molliou 31
Case postale
1131 Tolochenaz
Switzerland
Tel: +41 (0) 21 802 70 00
Fax: +41 (0) 21 802 79 00

UC201910267EE © 2019 Medtronic.
All Rights Reserved.

medtronic.eu